




Your supplier of high pressure laboratory instruments

AI-70 ABRASIVITY INDENTERS

ASTM D7625-10

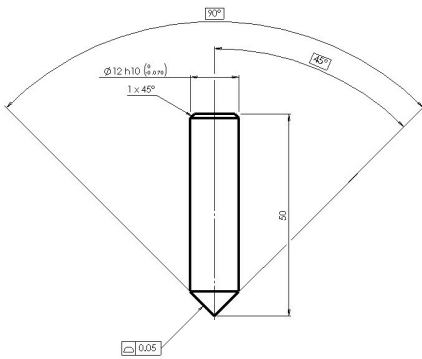

Description

The indenter serves as stylus pin for original and west Cerchar abrasivity index testers. The indenter is designed as per the ASTM D7625 method for determining the rock abrasiveness. The tip of test pin comprises a conical angle of 90°. Typical lifetime is around 10 tests per unit.

Specifications

Standard	ASTM D7625-10.
Diameter	12-mm
Length	50-mm
Cone	90°
Rockwell hardness	HRC 54/56 (200 DaN/mm ²)

Benefits

ASTM-compliant, original version

Fast delivery